



UNIVERSAL SCIENTIFIC

6210 Campbell Drive. Madison, OH 44057-2003
PHONE: 1-440-428-7800 or 1-440-829-8961 FAX: 1-440-428-8650
E-MAIL: tomheckman@universalscientificinc.com

Model 5434SF

Freestanding Spindle Rack Glassware Washer



Product Description

The 5434SF is perfect for washing and drying narrow-neck and general purpose glassware. Adjustable cycle programs give the user full control in setting water temperature, wash and rinse times, and number of rinses. Temperatures reaching up to 82° C (180° F) may be programmed to sanitize glassware. The result is thoroughly clean, dry labware, ready for use in sensitive research.

The washer is outfitted with a Premium upper basket with fold down shelves to wash your everyday bowls, beakers, and medium sized labware. Meanwhile, the bottom has 34 Spindles for washing narrow-neck labware. In addition, a premium lower basket is included and can be used instead of the Lower Spindle Rack.

Save Water - Save Electricity - Save Time with the 5434SF!

www.universalscientificinc.com

Product Dimensions

Height (Adjustable)	35 1/8" – 36 3/4"
Width	25 1/8"
Depth	28 7/8"
Depth with Door Open 90°	55 1/2"

Features

- Hidden Control Panel w/LCD Display
- Extra Tall Tub Interior Design
- 6 Temperatures Settings – Up to 82° C/180° F
- Quiet 48 dB(A) Operation
- 1 – 24 Hours Delayed Start
- End-of-Program Signal
- Stainless Steel Triple Filtration Self-Cleaning System
- Turbo Drying System
- Program Memory Recalls Last Program

Cycles

6 Cycles – Daily, Heavy, Normal, Quick, Eco, Rinse and Hold

Rack System

Adjustable Shelves

Premium Upper Basket – 9" Clearance

Premium Lower Spindle Rack – 34 customizable length spindles. 8" standard

Additional Lower Basket – Option to switch out Spindle Rack to wash larger wide mouth labware. 13 3/4" Clearance

Quality

- 18:9 Stainless Tank, Filter and Spray Arms
- Anti-Jam, Auto-Reversing Drain Pump
- Balanced Door
- Check Valve
- Galvanized Steel Base Pan
- PEX (cross-linked polyethylene) Water Line
- 6 foot long power cord included


Specifications

Stainless Steel Hidden Heating Element	
Adjustable leveling legs	
Annual Energy Consumption	231 kWh
Water Usage (Gallons)	.7 per fill
Voltage Rating	120V/60Hz/15A
Shipping Weight (Approx.)	195 Lbs.

Available Colors



White

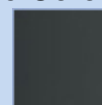

Black



Stainless